



The Glide N' Glow


Please read this instruction manual carefully before use.



1. Press and hold to turn on/off. Double-click to adjust the strength up 3 levels. Try using your thumbnails to double click for strength levels.
2. The device is set for five minutes. At four minutes you will hear a double beep letting you know that there is one minute left.
3. When the battery is low, the three green mode indicators automatically flash 10 times. The device will then power off automatically.
4. The device will not work while charging. The intensity indicator lights will flash while charging and go solid green when fully charged.


| | |
|--|--------------------------|
| Product Name: Neck Beauty Device | Rated Voltage: 3.7V |
| Model No.: DS-8822 | Material: ABS + Plating |
| Unit Size: 85*106*45MM | Rated Input: 1.5W |
| Charging Time: 3.5hrs | Battery Capacity: 500mAh |
| Main Functions: Lifting and Tightening | |

Blue light

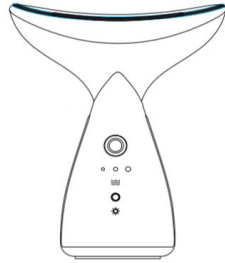


Mode 1: Vibration and blue LED. Use for calming inflammation.
(Blue LED 460-465 nanometers)

Purple light

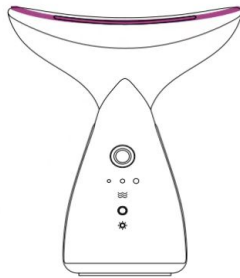

Mode 2: Vibration, EMS (activates collagen and elastin for toning) and purple LED. This is a combination of red and blue LED together. Improves oxygen to cells and calms the skin. Great for toning.

Red light

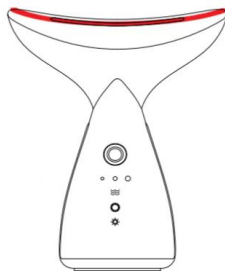

Mode 3: Vibration, heat (45 degrees) and red LED accelerates blood circulation, promotes growth of collagen and elastin.
(Red LED 620-625 nanometers)

Great for anti-aging.

Product Features

The streamlined head of this device is specially designed for the neck but works very well over the full face. See video.

There is sonic vibration with all three modes. This vibration stimulates blood circulation and metabolism, reduces edema, and creates relaxation. This is all effective in helping to promote anti-aging.

The three modes are: calming (blue light), toning (purple light), heat (red light),

There are three levels of intensity for each mode. Double click on/off button to raise intensity. Try using your thumbnails to double click.

How to Use

Clean face and rinse well. Apply a skin care product with plenty of moisture. Recommended products are QUL Anti-aging moisturizer with a few drops of GLO Serum or Hyaluronic Primer. Use the Hydrating Rose Water Spray to rehydrate if mixture gets dry.

We recommend applying products to one area at a time, working that area, before applying products to the next area. Example: Apply product mixture to the neck and work the neck for 5 minutes. Then apply to half of face and forehead. After doing 5 minutes, apply to the second half of the face.

Turn on the unit and choose the mode according to your skincare needs. Use slow upward movement.

Charging Info

Use the USB cord to fully charge the device before using it. This takes about 3.5 hours. When fully charged it can be used continuously for up to 12 hours. You will see flashing green lights when it is charging and solid green lights when charging is complete.

If the product is not going to be used for a long period of time, please charge it every three months. This will extend the life of the device.

Misc. Info

Wipe the device with an anti-bacterial wipe being careful not to get fluid inside of the device. Make sure that the device is completely clean of product after each use.

The device is not waterproof. Do not place it in water or use in shower.

Contraindications: pregnancy, pacemaker, open wounds.

Video: <https://youtu.be/rfhtsx1NY7w>